


DuraFlow 2


CHRONIC DIALYSIS CATHETER

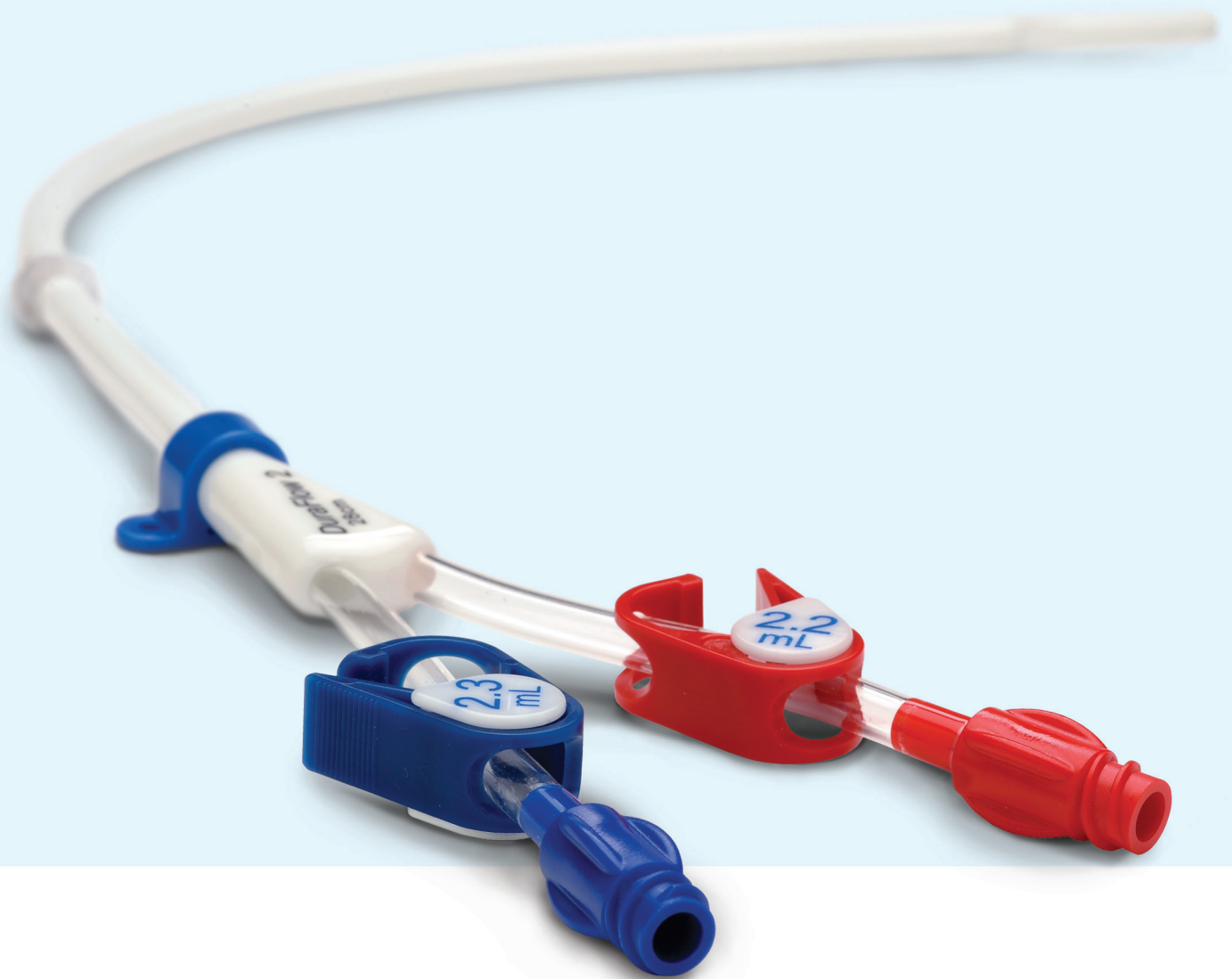

DuraFlow 2

CHRONIC DIALYSIS CATHETER

Engineered for Durability

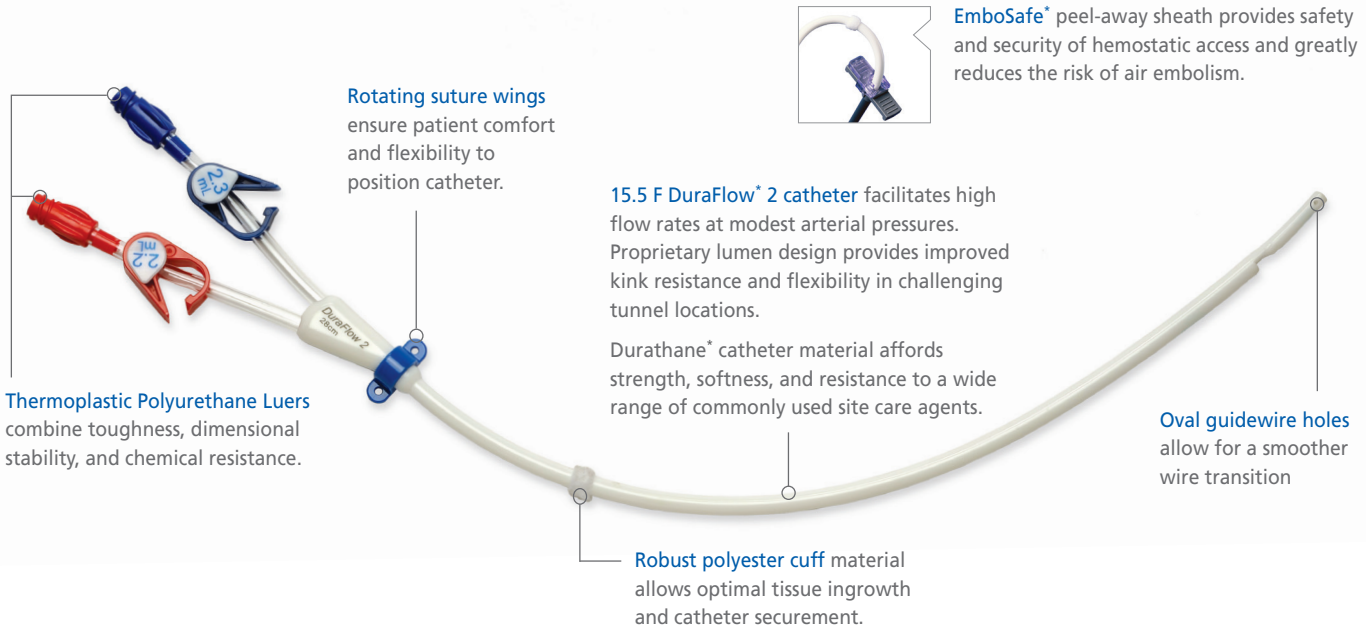

Rotating suture wings
ensure patient comfort
and flexibility to
position catheter.

Thermoplastic Polyurethane Luers
combine toughness, dimensional
stability, and chemical resistance.

15.5 F DuraFlow² catheter facilitates high
flow rates at modest arterial pressures.
Proprietary lumen design provides improved
kink resistance and flexibility in challenging
tunnel locations.

Durathane[®] catheter material affords
strength, softness, and resistance to a wide
range of commonly used site care agents.

Oval guidewire holes
allow for a smoother
wire transition

Robust polyester cuff material
allows optimal tissue ingrowth
and catheter securement.

EmboSafe[®] peel-away sheath provides safety
and security of hemostatic access and greatly
reduces the risk of air embolism.

Ordering Information (5 units per box)

| Description | Basic Kit | VascPak Kit | Catheter Only |
|-------------|--|--|--|
| 20 cm | <input type="checkbox"/> H787103012145 | — | — |
| 22 cm | <input type="checkbox"/> H787103012125 | — | — |
| 24 cm | <input type="checkbox"/> H787103012015 | <input type="checkbox"/> H787103012185 | <input type="checkbox"/> H787103017015 |
| 28 cm | <input type="checkbox"/> H787103012025 | <input type="checkbox"/> H787103012195 | <input type="checkbox"/> H787103017025 |
| 32 cm | <input type="checkbox"/> H787103012035 | <input type="checkbox"/> H787103012205 | <input type="checkbox"/> H787103017035 |
| 36 cm | <input type="checkbox"/> H787103012045 | <input type="checkbox"/> H787103012215 | <input type="checkbox"/> H787103017045 |
| 40 cm | — | — | <input type="checkbox"/> H787103017055 |
| 48 cm | — | — | <input type="checkbox"/> H787103017155 |
| 55 cm | — | — | <input type="checkbox"/> H787103017145 |

Ordering Information (1 unit per box)

| Description | Basic Kit | VascPak Kit |
|-------------|--|--|
| 24 cm | <input type="checkbox"/> H787103012011 | — |
| 28 cm | <input type="checkbox"/> H787103012021 | — |
| 32 cm | <input type="checkbox"/> H787103012031 | — |
| 36 cm | <input type="checkbox"/> H787103012041 | — |
| 40 cm | <input type="checkbox"/> H787103012055 | <input type="checkbox"/> H787103012225 |
| 48 cm | <input type="checkbox"/> H787103012155 | <input type="checkbox"/> H787103012235 |
| 55 cm | <input type="checkbox"/> H787103012115 | <input type="checkbox"/> H787103012245 |

BASIC KIT INCLUDES: (1) DuraFlow 2 Dialysis Catheter, (1) 16F Valved Breakaway Introducer, (1) 12F Dilator, (1) 14F Dilator, (1) Tri-Ball Tunneler with Sleeve, (2) Injection Caps, (1) #11 Blade Safety Scalpel, (1) 0.038" J/Flex Guidewire, (1) Adhesive Wound Dressing, (1) 18 Ga. x 2-3/4" Introducer Needle

VASCPAK[®] KIT INCLUDES: (1) DuraFlow 2 Dialysis Catheter, (1) Adhesive Wound Dressing, (1) 16F Valved Peelable Sheath, (1) Tri-ball Tunneler with Sleeve, (2) Injection Caps, (1) 12F Vessel Dilator, (1) 14F Vessel Dilator

IMPORTANT RISK INFORMATION

INDICATIONS FOR USE: The DuraFlow 2 Hemodialysis Catheter is indicated for use in attaining long-term vascular access for hemodialysis and apheresis. It may be inserted percutaneously and is primarily placed in the internal jugular vein of an adult patient. Alternate insertion sites include subclavian vein as required. Catheters greater than 40 cm are intended for femoral vein insertion. This catheter is indicated for > 30 days long term placement. The valved peelable introducer is

intended for use in the percutaneous insertion of catheters in the venous system.

CONTRAINDICATIONS: The catheter is intended for long-term vascular access only and should not be used for any purpose other than indicated in these instructions. The valved peelable introducer sheath is not designed for use in the arterial system.

WARNINGS AND PRECAUTIONS: Please see package insert for complete list of warnings and precautions

POTENTIAL COMPLICATIONS: Consult package insert for a complete list of potential complications.

CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician.

Indications, contraindications, warnings and instructions for use can be found in the instructions for use supplied with each device. Observe all instructions prior to use. Failure to do so may result in patient complications.



USA > 14 Plaza Drive, Latham, NY 12110 > tel: 800-772-6446 or 518-798-1215 > fax: 518-798-1360
International > Haaksbergweg 75 (Margrietoren), 1101 BR, Amsterdam Z-O > The Netherlands
tel: +31 (0)20 753 2949 > fax: +31 (0)20 753 2939

www.angiodynamics.com

*AngioDynamics, the AngioDynamics logo, DuraFlow, EmboSafe, Durathane and Vascpak are trademarks and/or registered trademarks of AngioDynamics, Inc., an affiliate or a subsidiary.

© 2015 AngioDynamics, Inc. MLC 456 US Rev E 04/15