

# MINI STICK MAX

## Coaxial Microintroducer Kit

**STARTING HAS NEVER BEEN EASIER.**

AngioDynamics is committed to providing physicians the best tools possible for the best patient outcomes, no matter how big or small the device. Designed to start each procedure with ease and efficiency, the AngioDynamics Mini Stick MAX Coaxial Microintroducer Kits incorporate proprietary designs and features that aid in simplifying set up and gaining vascular access.

### NEEDLE OPTIONS

4 cm and 7 cm length  
available with or without  
echogenic tip

### SHEATH DILATOR OPTIONS

10 cm and 15 cm  
sheath dilator options

### CONTAINMENT CLIP

designed to keep all components  
together in a sterile field

### THREE DIFFERENT 0.018" WIRE CONFIGURATIONS

- Stainless/Tungsten
- Nitinol/Tungsten
- Stainless Steel/Stainless Steel

**TRANSLUCENT HUB**  
enables visualization  
of flash back

### INTERLOCKING COAXIAL DESIGN

in 4F and 5F with stiff  
or standard dilator




angiodynamics

GL/VA/BR/1973 Rev 01 08/2023

**TAPERED DILATOR**  
for a smooth transition at the tip

**PROPRIETARY COAXIAL DESIGN**  
for enhanced ease of insertion

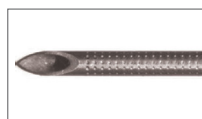

**DIMPLING ON DISTAL TIP**  
for enhanced echogenicity

**BLACK BEVEL INDICATOR**  
provides easy indication of  
bevel orientation

KIT COMPONENTS							MINI STICK MAX COAXIAL MICROINTRODUCER KITS <sup>1,2</sup>				BOX QTY
0.018" WIRE		21G NEEDLE				SHEATH	STANDARD		STIFFENED		
CORE/TIP MATERIAL	LENGTH (CM)	ECHOGENIC TIP	LENGTH (CM)	BEVEL	HUB	LENGTH (CM)	4F	5F	4F	5F	
NITINOL / TUNGSTEN	45	YES	4	B	GREEN	10	H965457711	H965457751	H965457811	H965457851	10
			7	B	GREEN	10	H965457501	H965457561	H965457531	H965457591	10
		NO	4	B	GREEN	10	H965457691	H965457731	H965457781	H965457831	10
				B	GREEN	15				H965457861	10
			7	B	GREEN	10	H965457631	H965457671	H965457621	H965457661	10
				B	GREEN	15	-	-	-	H965457761	10
STAINLESS STEEL / STAINLESS STEEL	45	YES	4	B	GREEN	10	H965457701	H965457741	H965457791	H965457841	10
			7	B	GREEN	10	H965457481	H965457541	H965457511	H965457571	10
		NO	4	B	GREEN	10	H965457681	H965457721	H965457771	H965457821	10
			7	B	GREEN	10	H965457611	H965457651	H965457601	H965457641	10
STAINLESS STEEL/ TUNGSTEN	45	YES	7	B	GREEN	10	H965457491	H965457551	H965457521	H965457581	10

Each kit includes: one .018" guidewire, one 10 cm sheath dilator/introducer and one 21 gauge needle. † 15 cm sheath

#### FOOTNOTES:

<sup>1</sup>Not all products may be available for sale in your market. Please consult with your local AngioDynamics Partner for product availability and additional product information.

<sup>2</sup>Low volume configurations may be built to order and require a greater lead time.

#### IMPORTANT RISK INFORMATION

##### MINI STICK MAX COAXIAL MICROINTRODUCER KIT

**INDICATION FOR USE:** The coaxial microintroducer kit is used for the percutaneous introduction of a guidewire into the vascular system.

**CONTRAINDICATIONS:** None known.

Refer to Directions for Use and/or User Manual provided with the product for complete Instructions, Warnings, Precautions, Possible Adverse Effects and Contraindications prior to use of the product.

**CAUTION:** Federal (USA) law restricts this device to sale by or on the order of a physician.

USA > 14 Plaza Drive, Latham, NY 12110 > Tel: 800-772-6446 or 518-798-1215 > Fax: 518-798-1360  
International > Haaksbergweg 75 (Margrietoren), 1101 BR, Amsterdam Z-O > The Netherlands  
Tel: +31 (0)20 753 2949 > Fax: +31 (0)20 753 2939  
www.angiodynamics.com



\*AngioDynamics, the AngioDynamics logo and Mini Stick are trademarks and/or registered trademarks of AngioDynamics, Inc., an affiliate or a subsidiary. All other trademarks are property of their respective owners. © 2023 AngioDynamics, Inc. GL/VA/BR/1973 Rev 01 08/2023